

Care for A Better Life

Suzhou Kangyouwei Medical Technology Co., Ltd.

Web: www.zodelo.com

TEL: +86-13776135860

Post Code: 215000

Email: info@zodelo.com

Address: Building 9, No. 28 Xingzhong Road, Wuzhong
Economic Development Zone, Suzhou, China



Zodelo Ostomy Care Products



CONTENTS

P 02 • Company Introduction

R&D Investment	P 03
Company Honors	P 04
Company Culture	P 04

P 05 • Zodelo Products

Airlora Series	P 06
Calmlora Series	P 08
Mellora Series	P 17
Joylora Series	P 25
Stoma Accessories	P 27

P 31 • Customization Process

Product analysis	P 31
Customization Process	P 32

P 33 • Product Structure

Sealing Method	P 33
Scale Substrate	P 33
EVOH or PVDC film	P 34
Fabric/back window design	P 34
Exhaust valve	P 35
Carbon filter	P 35
Medical hydrolloid	P 36
Spunlace extender	P 36
Bag Shape	P 37
Two-piece system with flange	P 37
Anti-backflow design	P 38
Logo customization	P 38

COMPANY INTRODUCTION

Vison:

Care for a better life.



About US:

Suzhou Kangyouwei Medical Technology Co., Ltd. was founded in 2007. It is an innovation-driven medical company focused on the development, manufacturing, and sales of ostomy and related medical products. The company is dedicated to building its own brand, ZODELO.

The company is located in Suzhou, China. Its manufacturing facility covers approximately 10,000 square meters. Kangyouwei employs over 200 skilled technicians and more than 40 R&D engineers. The annual production capacity reaches 65 million units, with an annual output value of around RMB 150 million.

From 2021 to 2025, the ZODELO ostomy pouch series ranked No. 1 in sales on major e-commerce platforms in China. This reflects the brand's reliable quality, strong market acceptance, and growing user trust.

In June 2025, ZODELO began its global expansion and started building an international distributor network. Helping more users live better lives is the long-term mission of ZODELO.

R&D Investment



Kangyouwei Medical Technology Co., Ltd. regards innovation as the lifeblood of its survival and development, continuously focusing on user needs to drive product development and upgrades.

As of 2025, our company has been granted 50 patents, including 3 invention patents and 46 utility model patents.

Kangyouwei is ISO 13485 certified and has completed FDA registration and CE self-declaration in accordance with applicable regulations.

Company Honors



Title of Private
Technology Enterprise



National High-Tech
Enterprise Certification



Enterprise Credit
Management System
Certification


Company Culture

- Our Slogan: Care for a better life
- Our Focus: We focus not only on products, but on every user's real experience.
- Our Quality: We believe quality is not tested — it is built. Every product we make is a tribute to life.

Zodelo Five Major Product Series

Airlora Series

Entry-Level Model

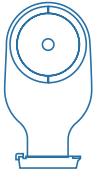

Calmlora Series

Basic Model



Mellora Series

Upgraded Model



Joylora Series

Premium Model



Stoma Accessories

Airlora Series

— Entry-Level Model



Airlora Series

■ One-Piece Ostomy Bag



1.

Double-layer baseplate
Strong adhesion
2.

Filter odors
Prevent ballooning
3.

Double-Sided Fabric
Comfortable wear
4.

Locks in Odors
Worry-FreeSocializing
5.

Hydrocolloid Baseplate
Soft and Gentle
6.

Precise Scale
Easy to Cut

Model	214047
Max cut	65mm
Skin barrier diameter	130mm
Fabric	Double sides brown non-woven fabric
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Bag color and Capacity	brown,350ML
Closure	Clamp
Other structures	Front window,small arc-shaped carbon piece
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box

Calmlora Series

——Basic Model



Calmlora Series

■ One-Piece Flat Ostomy Bag



1. Integrated vent & clean design
2. Reinforced ring twice as secure
3. Double-hook Velcro for enhanced durability
4. Easy-pull outlet design
5. Strong seal locks in odors
6. Double-sided non-woven fabric for greater comfort

Model	214512
Max cut	65mm
Skin barrier diameter	130mm
Fabric	Double sides white non-woven fabric
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Bag color and Capacity	Double sides White, 350ML
Closure	Velcro
Other structures	Front window, blue exhaust valve
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box, 10 boxes/ctn

Calmlora Series

■ One-Piece Convex Ostomy Bag



1. Integrated vent & clean design
2. Reinforced ring twice as secure
3. Double-hook Velcro for enhanced durability
4. Easy-pull outlet design
5. Strong seal locks in odors
6. Double-sided non-woven fabric for greater comfort

Model	214518
Max cut	45mm
Skin barrier diameter	130mm
Fabric	Double sides white non-woven fabric
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Bag color and Capacity	Double sides White, 350ML
Closure	Velcro
Other structures	Front window, blue exhaust valve
Box size	
Carton size	370x250x280mm
Carton specification	10pcs/box, 10 boxes/ctn

Calmlora Series

■ Two-Pieces Ostomy Bag

-Transparent bag



1. Two-Piece Design
Easy Use & Cleaning

2. Integrated vent &
clean design

3. Soft Buckle for
Comfortable Wear
4. Compatible with Two
Baseplates for Different
Ostomy Needs

5. Easy-pull outlet design

6. Double-hook Velcro for
enhanced durability

Model	214531
Flange	60mm
Fabric	Single side white non-woven fabric
Bag color	Single side white
Other structures	Blue exhaust valve
Closure	Velcro
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

■ Two-Pieces Ostomy Bag



1. Two-Piece Design
Easy Use & Cleaning

2. Integrated vent &
clean design

3. Soft Buckle for
Comfortable Wear
4. Compatible with Two
Baseplates for Different
Ostomy Needs


5. Easy-pull outlet design

6. Double-hook Velcro for
enhanced durability

Model	214541
Flange	60mm
Fabric	Double sides white non-woven fabric
Bag color	Double sides white
Other structures	Front window, blue exhaust valve
Closure	Velcro
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

■ One-Piece Disposable Closed-end Ostomy Bag

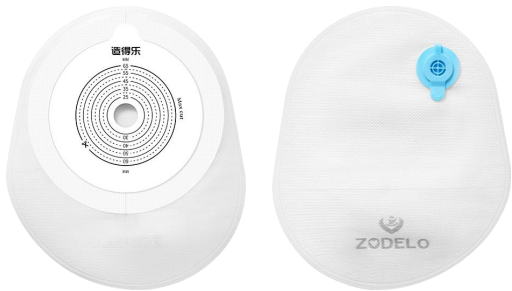

1. Integrated vent & clean design

2. Reinforced ring twice as secure

3. Strong seal locks in odors
4. Double-sided non-woven fabric for greater comfort

5. Precise Scale Easy to Cut

6. Double-layer baseplate Strong adhesion

Model	214812
Max cut	65mm
Skin barrier diameter	130mm
Fabric	Double sides white non-woven fabric
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Bag color and Capacity	Double sides White,200ML
Closure	Closed-end
Other structures	Front window, blue exhaust valve
Box size	
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

■ Two-Piece Disposable Closed-end Ostomy Bag



1. Two-Piece Design Easy Use & Cleaning

2. Integrated vent & clean design

3. Soft Buckle for Comfortable Wear
4. Compatible with Two Baseplates for Different Ostomy Needs

5. Precise Scale Easy to Cut

6. Double-sided non-woven fabric for greater comfort

Model	215021
Flange	55mm
Fabric	Double sides white non-woven fabric
Bag color	Double sides white
Other structures	Front window, blue exhaust valve
Closure	Velcro
Box size	
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

■ One-Piece Urostomy Bag



1.

Hypoallergenic High-Adhesive Hydrocolloid
2.

Anti-Reflux Infection Preventio
3.

Clear Scale for Easy Cutting
4.

Upgraded Baseplate for Leak Prevention
5.

Breathable Non-Woven Lining Skin-Friendly
6.

Soft Plug Seal, Leak-Proof Design

Model	113042
Max cut	65mm
Skin barrier diameter	130mm
Fabric	Single side white non-woven fabric
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Bag color	Single side transparent
Closure	Drain valve
Other structures	Anti-backflow design
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

■ Two-Pieces Urostomy Bag



1.

Two-Piece Design Easy Use
2.

Anti-Reflux Infection Preventio
3.

Soft Buckle for Comfortable Wear
4.

Compatible with Two Baseplates for Different Ostomy Needs
5.


Breathable Non-Woven Lining Skin-Friendly
6.

Soft Plug Seal, Leak-Proof Design

Model	113412
Flange	60mm
Fabric	Single side white non-woven fabric
Bag color	Single side transparent
Other structures	Anti-backflow design
Closure	Drain valve
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Calmlora Series

Two-Pieces Flat Skin Barrier

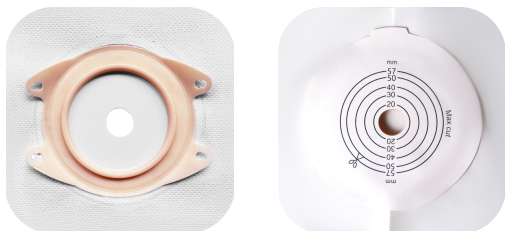

1. Square spunlace reinforcement
2. Precise Scale Easy to Cut
3. Double-layer baseplate Strong adhesion
4. Single-piece PVC blister packaging
5. Enhanced hydrocolloid ideal for sweaty skin
6. Hydrocolloid Baseplate Soft and Gentle

Model	02833
Max cut	57mm
Skin barrier size	130mm
Flange	60mm
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Printing line base material	Release paper
Box size	170x160x40mm
Carton size	
Carton specification	5 pcs/box,

Calmlora Series

Two-Piece Convex Skin Barrier



1. Outer spunlace adhesive sweat-resistant, breathable for a stronger hold
2. Precise Scale Easy to Cut
3. Double-layer baseplate Strong adhesion
4. Single-piece PVC blister packaging
5. Enhanced hydrocolloid ideal for sweaty skin
6. Hydrocolloid Baseplate Soft and Gentle

Model	11037
Max cut	45mm
Skin barrier size	130mm
Flange	60mm
Barrier material	Medical hydrocolloid and Non-woven reinforced adhesive ring
Printing line base material	Release film
Box size	170x160x70mm
Carton size	
Carton specification	5 pcs/box,

Mellora Series

— Upgraded Model



Mellora Series

■ One-Piece Ostomy Bag



- 1. Softer baseplate for better fit and no edge lifting
- 2. Upgraded Adhesion Longer Lasting Use
- 3. Oval Adhesive Baseplate Strong Absorbency & Durability
- 4. Chinese-style design elegant and comfortable
- 5. Upgraded Gray Film for Enhanced Odor Barrier
- 6. Integrated Venting & Cleaning Design Convenient and Practical

Model	214513
Max cut	70mm
Skin barrier diameter	120*130mm
Fabric	Double sides grey non-woven fabric
Barrier material	Medical hydrocolloid
Bag color	Grey
Closure	Velcro
Other structures	Front window,Grey Exhaust Valve
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Mellora Series

Two-Pieces Ostomy Bag

-Transparent bag



透明二件式袋
子无法喷印
LOGO

- 1. Two-Piece Design
Easy Use & Cleaning
- 2. Integrated vent &
clean design
- 3. Reinforced Edging
Long-Lasting Adhesion
- 4. Compatible with Two
Baseplates for Different
Ostomy Needs
- 5. Secure Lock
Easy On/Off
- 6. High-Barrier PVDC for
Superior Odor Isolation

Model	214563
Flange	60mm
Fabric	Single side brown non-woven fabric
Bag color	Brown and transparent
Other structures	Brown Exhaust Valve
Closure	Velcro
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box, 10 boxes/ctn

Mellora Series

Two-Pieces Ostomy Bag



- 1. Two-Piece Design
Easy Use & Cleaning
- 2. Integrated vent &
clean design
- 3. Reinforced Edging
Long-Lasting Adhesion
- 4. Compatible with Two
Baseplates for Different
Ostomy Needs
- 5. Secure Lock
Easy On/Off
- 6. High-Barrier PVDC for
Superior Odor Isolation

Model	214564
Flange	60mm
Fabric	Double sides brown non-woven fabric
Bag color	Double sides brown
Other structures	Brown Exhaust Valve
Closure	Velcro
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box, 10 boxes/ctn

Mellora Series

■ One-Piece Disposable Closed-end Ostomy Bag



1. Softer baseplate for better fit and no edge lifting

2. Upgraded Adhesion Longer Lasting Use

3. Oval Adhesive Baseplate Strong Absorbency & Durability
4. Chinese-style design elegant and comfortable

5. Upgraded Gray Film for Enhanced Odor Barrier

6. Integrated Venting & Cleaning Design Convenient and Practical

Model	214813
Max cut	70mm
Skin barrier diameter	132mm
Fabric	Double sides grey non-woven fabric
Barrier material	Medical hydrocolloid
Bag color	Grey
Closure	Closed-end
Other structures	Front window, grey exhaust valve,200ML
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Mellora Series

■ Two-Piece Disposable Closed-end Ostomy Bag



1. Two-Piece Design Easy Use & Cleaning

2. Integrated vent & clean design

3. Soft Buckle for Comfortable Wear
4. Compatible with Two Baseplates for Different Ostomy Needs


5. Upgraded Gray Film for Enhanced Odor Barrier

6. Double-sided non-woven fabric for greater comfort

Model	214571
Flange	60mm
Fabric	Double sides grey non-woven fabric
Bag color	Grey
Other structures	Front window, grey exhaust valve,200ML
Closure	Closed-end
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box,10 boxes/ctn

Mellora Series

Two-Pieces Flat Skin Barrier

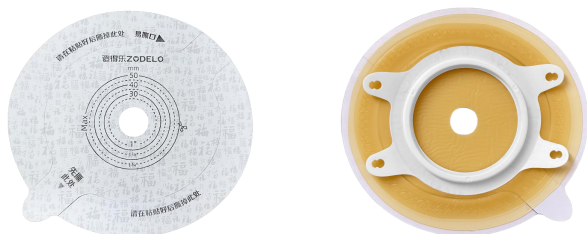

1. Upgraded Self-Development Adhesive
2. Precise Scale Easy to Cut
3. Double-layer baseplate Strong adhesion
4. Single-piece PVC blister packaging
5. Sweat-Resistant:With Sweat -Proof Function
6. Hydrocolloid Baseplate Soft and Gentle

Model	60132
Max cut	50mm
Skin barrier size	120*130mm
Flange	60mm
Can be manufactured	
Barrier material	Medical hydrocolloid
Printing line base material	Release paper
Box size	
Carton size	
Carton specification	5pcs/box,

Mellora Series

Two-Pieces Convex Skin Barrier



1. Upgraded Self-Development Adhesive
2. Precise Scale Easy to Cut
3. Double-layer baseplate Strong adhesion
4. single-piece PVC blister packaging
5. Sweat-Resistant:With Sweat -Proof Function
6. Hydrocolloid Baseplate Soft and Gentle

Model	11036
Max cut	45mm
Skin barrier size	120*130mm
Flange	60mm
Can be manufactured	
Barrier material	Medical hydrocolloid
Printing line base material	Release paper
Box size	
Carton size	
Carton specification	5pcs/box,



Joylora Series

— Premium Model

Joylora Series

Two-Pieces Ostomy Bag



- 1. Softer baseplate for better fit and no edge lifting
- 2. Upgraded Adhesion, Longer Lasting Use
- 3. Oval Adhesive Baseplate Strong Absorbency & Durability
- 4. Chinese-style design elegant and comfortable
- 5. Upgraded Gray Film for Enhanced Odor Barrier
- 6. Integrated Venting & Cleaning Design Convenient and Practical

Model	214562
Flange	60mm
Fabric	Double sides yellow non-woven fabric
Barrier material	Medical hydrocolloid
Bag color	yellow
Closure	Velcro
Other structures	Front window, large crescent-shaped carbon piece
Box size	220x170x50mm
Carton size	370x250x280mm
Carton specification	10pcs/box, 10 boxes/ctn

Stoma Accessories



Stoma Accessories



1. Ostomy powder

30g



2. Skin barrier film

30ml 60ml 100ml



3. Skin barrier film (sheet wipe)

30pcs/box



4. Ostomy paste

60g



5. Skin barrier ring

2mm 4mm 10pcs/box



6. Adhesive remover

30ml 60ml 100ml

Stoma Accessories



7. Elastic barrier strips

spunlace fabric


hydrocolloid



8. Adhesive remover
(sheet wipe)

20pcs/box

Stoma Accessories

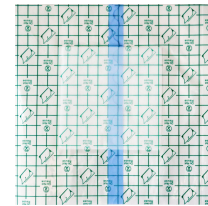

13. Waterproof Ostomy
Shower Film



9. Ostomy support belt

regular

summer lightweight



10. Ostomy belt




11. Stoma curved scissors



12. 500ml washing bottle

ONE-PIECE PRODUCT ANALYSIS PROFILE

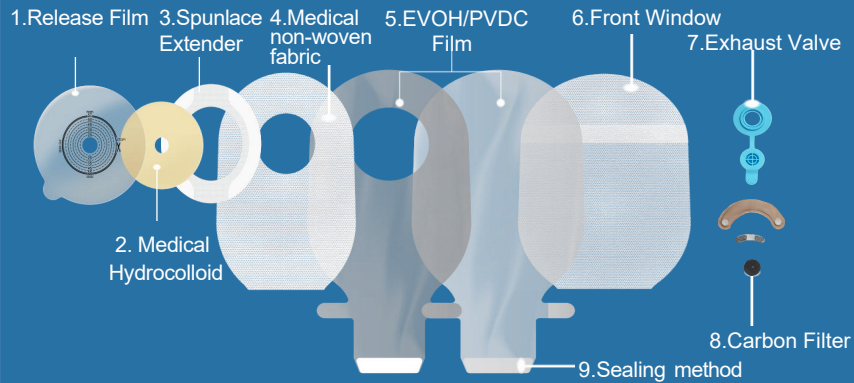

1. **Release film**(Strong adhesive, easy to peel off)
2. **Medical hydrocolloid**(Good stickiness, long-lasting and skin-friendly, no residue on the skin)
3. **Spunlace extender**(Reinforce the baseplate, prevent warping, not prone to allergies)
4. **Medical non-woven fabric**(Soft, comfortable and ventilate, reduce friction noise, avoid embarrassment)
5. **EVOH film**(7-layer film co-extruded, strong and durable)
6. **Back window**(Double layer medical non-woven fabric, reduce friction sound, more concealed)
7. **Exhaust valve**(Quickly discharge the flatulence in the bag)
8. **Carbon filter**(Good adsorption, effectively filter excrement odor)
9. **Sealing method**(The sealing operation is convenient and the sealing is tight)




Ostomy Bag Customization Process

If you want to have a unique ostomy bag, you can customize it as you like.

You can choose step by step according to the following process, and the corresponding structure is introduced later.

Please choose 1~8 for one piece system, please choose 1~9 for two-piece system.

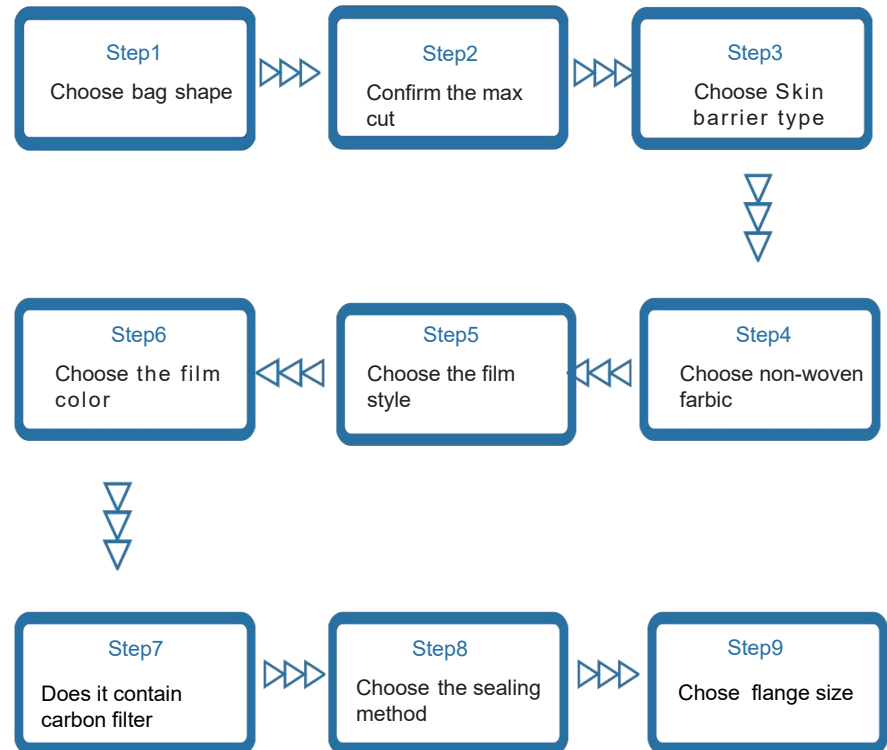

Structure 1



Sealing Method

4 types of sealing methods:

3 types for **colostomy bag**, namely---Clamp, Hook and Loop, Velcro and Medical Twist Tie

1 type for **urostomy bag**: Drain Valve.



Clamp



Velcro



Medical Twist Tie



Drain Valve

Structure 2



Scale Substrate

2 types of base material for cutting hole printing line:

Release Paper and Release Film, which have the same function, protecting the skin barrier and printing scale marks, and making it convenient for patients to cut.



Structure 3



EVOH or PVDC film

The body of the ostomy bag can be made of EVOH film or PVDC film.

EVOH film is lightweight, durable, and offers excellent gas barrier properties.

PVDC film excels in moisture resistance and odor control.



Structure 4



Fabric/back window design

3 types of fabrics used for ostomy bags

white non-woven fabric


gray non-woven fabric

brown non-woven fabric

3 types of fabric



window design

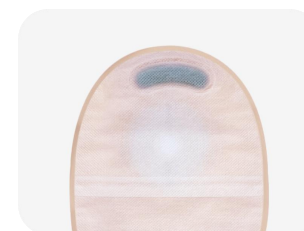

Structure 5



Exhaust valve

When the bag is full of gas, the gas release valve allows for quick and easy venting with simple operation.




Structure 6



Carbon Filter

Carbon pieces come in three different shapes: small black discs, small arcs, and large crescents. It absorbs odors and can also release gas, though venting is slower and requires punching holes.

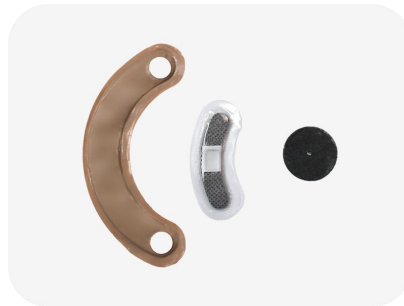

Structure 7



Medical hydrocolloid

The skin barrier is made of medical-grade hydrocolloid. It is translucent, sticky, gentle on the skin, comfortable to wear, and can last 5–7 days with proper use.



Structure 8



Spunlace extender

It enhances skin barrier adhesion, prevents warping, and keeps the bag secure—ideal for active patients.




Structure 9



Bag Shape

The bag comes in various shapes and supports **customization** to meet different needs.

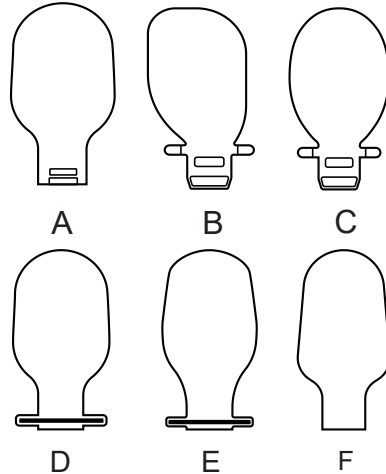

Structure 10



Two-piece system with flange

The flange system is made of translucent polyethylene and used with a buckle. It has 3 belt tabs for easy reinforcement and comes in three sizes: 45, 57, and 70.



Structure 11



Anti-backflow design

Used in urostomy bags, it prevents liquid backflow and protects the skin around the stoma.



Structure 12



Logo customization

We offer personalized customization options for the Ostomy baseplate, non-woven bag, and product packaging.

